



The Tradewind 2000 Alternating Pressure Mattress Replacement System

Active pressure relief for those Patients assessed at medium to high risk of developing pressure sores, or with superficial established tissue damage.

The Tradewind 2000 mattress replacement system is a value for money, general purpose dynamic support surface, providing effective pressure relief with durable and optimum Patient comfort. A true alternating support surface, providing effective pressure relief by gently inflating and deflating the individual cells alternately over a ten minute cycle.

Compact Unobtrusive Pump Design:





- Convenient moulded retractable carrying handle.
- Convenient moulded retractable hanging hooks designed to fit most beds.
- Ergonomic injection moulded enclosure design.
- Simple controls, manual comfort control setting with convenient patient weight guide for optimum Patient comfort and support.

Advanced operator features:

- Pressure control: Adjustable to Patient weight / Adjustable to Patient position / Maximum Patient weight 160Kg (352lbs) (25 Stones).
- Sophisticated Alarm Indicators: Audible and visual power fail alarm / Visual low pressure alert indicator / Integral battery backup for alarms.
- Transport facility: this facility is easy to activate by simply interconnecting the mattress snap lock couplings together.
- Lightweight and compact.
- Minimal noise and vibration levels.
- Tough durable casing.
- Eliminates recessed areas where dirt may accumulate.
- Illuminated on/off switch.
- Snap fit tube couplings for secure connection and to prevent accidental disconnection of the mattress tubing from the pump.

2000 Alternating Pressure Mattress System

- Fitted with quick release couplings, which can be interconnected, to equalise the mattress pressure, and keep the mattress inflated during a power failure or when the pump is not required.
- The Modular design and construction of the Tradewind 2000 ensures all cells can be individually replaced making it an extremely cost effective option.
- 17 individually replaceable PU (polyurethane) air cell construction for long lasting durability.
- Conforms to the contours of the body ensuring effective Patient support and comfort.
- Fitted with a quick release CPR facility for fast deflation of the entire mattress. The CPR can be activated with a single handed release operation deflating the torso section in less than 20 seconds.
 - Fully detachable cover for easy cleaning.
 - The cover is vapour permeable, ultra two way stretch, fire retardant, soft but robust and is water resistant.
 - Compact for storage.

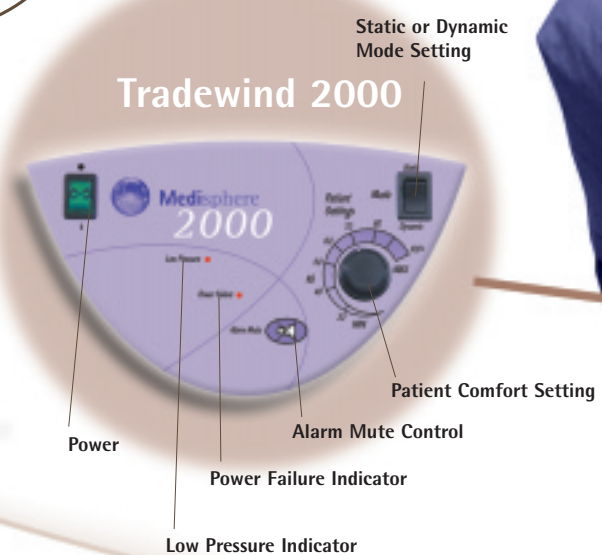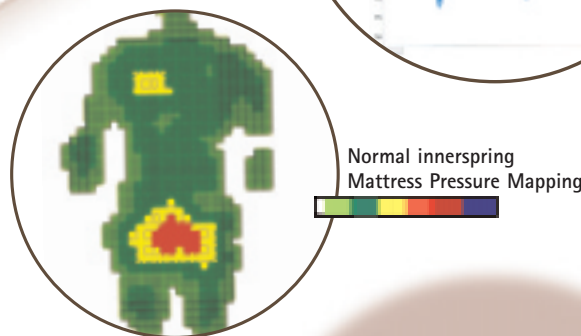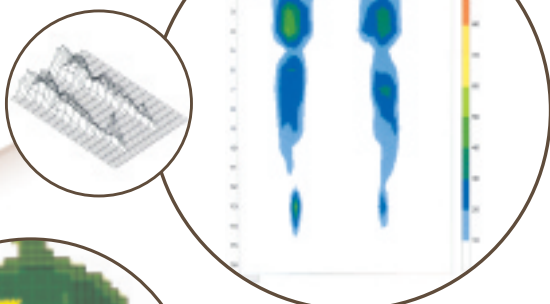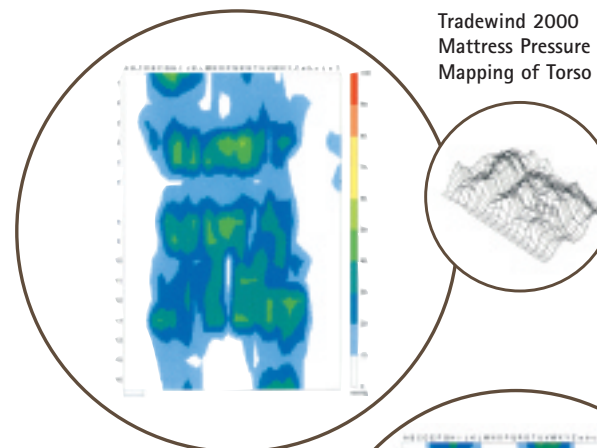

Fitted with quick release couplings, which can be interconnected, to equalise the mattress pressure, and keep the mattress inflated during a power failure or when the pump is not required.



The mattress cover is a multi-stretch material that is moisture vapour permeable, helping to reduce shear friction and maceration. The Tradewind 2000 Mattress has a removable bacteriostatic top cover that is easily cleaned between Patients or removed completely for laundering.



20 individually replaceable PU (polyurethane) air cell construction for long lasting durability.

Fitted with a quick release CPR facility for fast deflation of the entire mattress. The CPR can be activated with a single handed release operation deflating the torso section in less than 20 seconds.



Security. The mattress is fitted with adjustable straps to ensure secure fittings, to the bed frame.



Cost Effectiveness. The Modular design and construction of the Tradewind 2000 ensures all cells can be individually replaced making it an extremely cost effective option.



Compact For Storage. As storage space is at hospitals or equipment loan stores is at a premium the Tradewind 2000 is designed to take up minimum space when the system is not in use.



The therapeutic support surface is lightweight and conveniently rolls into its own carrying case.